


# Zematra Total Sediment Tester

ASTM D4870; IP375; ISO 10307

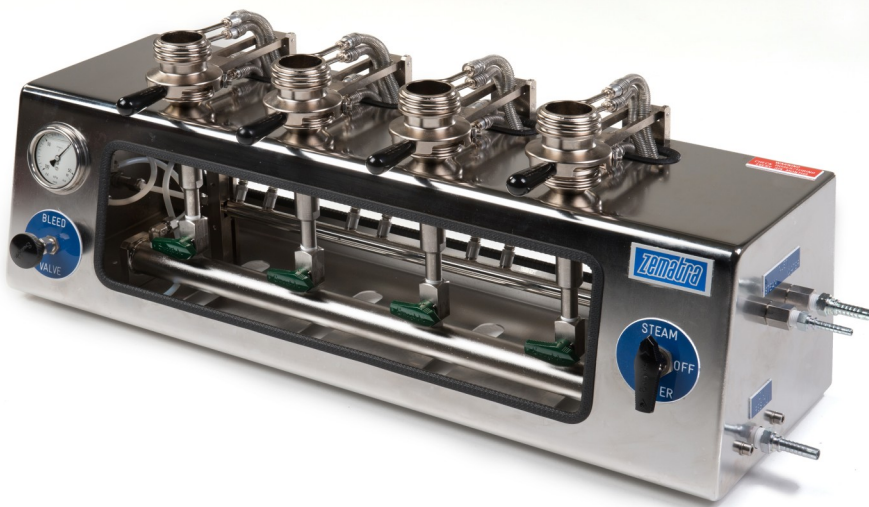

The Zematra 4 position Sediment Tester determines the insoluble material content of distillate and residual fuel oils.

It is a bench model which consists of four filtration cells, vacuum circuit with gauge and integrated cooling/heating circuit.

#### Optional:

- Vacuum pump
- Steam generator
- Waste container
- Thermometer
- Balance

#### Specifications:

- Measures sediment up to 0.5% m/m
- Samples viscosities up to 55 mm<sup>2</sup>/s at 100°C
- 10 +/- 0.5 gram sample size
- 4 test stations
- Dimension 72x30x30cm

Each cup is equipped with a sintered brass disc to support the microfibre filters.



Zematra, your reliable partner in testing



2139 High Tech Road | State College, PA | 16803 | USA  
800-676-6232 | 814-353-8000 | Fax: 814-353-8007  
sales@cannoninstrument.com | cannoninstrument.com